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Moist heat is a recognised, pleasant and well-tried method to treat and ease tension, pains and various complaints.

If you have any medical queries, please contact your physician or therapist.

Please read the instructions and safety notes prior to using the HYDROTHERM – heat pad:

- Do not use when going to sleep.
- Do not use on infants.
- Do not use if skin numb or poorly supplied with blood.
- Never use without the provided cover.
- Keep away from moisture.
- Individual treatments should not be longer than the preset 30 min.
- Do not pull cables or plugs.
- Disconnect from main power when not in use.
- Avoid damage to the control unit.
- Keep away from pointed objects.
- Do not use or put load on when folded.
- When used several times in succession, the pad cannot generate sufficient moisture. The special effect of moist heat is therefore not adequate. Simply allow the pad to cool before switching it back on. Afterwards the pad will function as usual.

Instructions for Cleaning and Maintenance

- Machine wash cover at 30°C
- Wipe pad carefully with a wet cloth. If necessary, add mild detergent.

Operation (see illustration on the right)

1. Switch on unit. The display performs an automatic check. The pad is in operation with the default settings.
2. Display: Settings control
3. For changing time and temperature, press selector switch and adjust settings with the +/- button.
4. The unit will switch off after the set time. Press OFF button to switch off manually (as in 1.).


Display / Indicator / Symbols

- Heater: Pad is heating up
- Clock: Time setting active
- Thermometer: Temperature setting active
- Key: + / - Key lock active
- Error: Indicator above error. Do not use pillow.
- Temperature: Temperature display in Celsius
- Time: Time display in minutes

Specifications

- Power consumption: 220 - 240 V / 50 Hz
- Temperature range: 30° - 85° C (65° preset)
- Duration: 1 - 60 min. (30 min. preset)

Operating Elements Control Unit

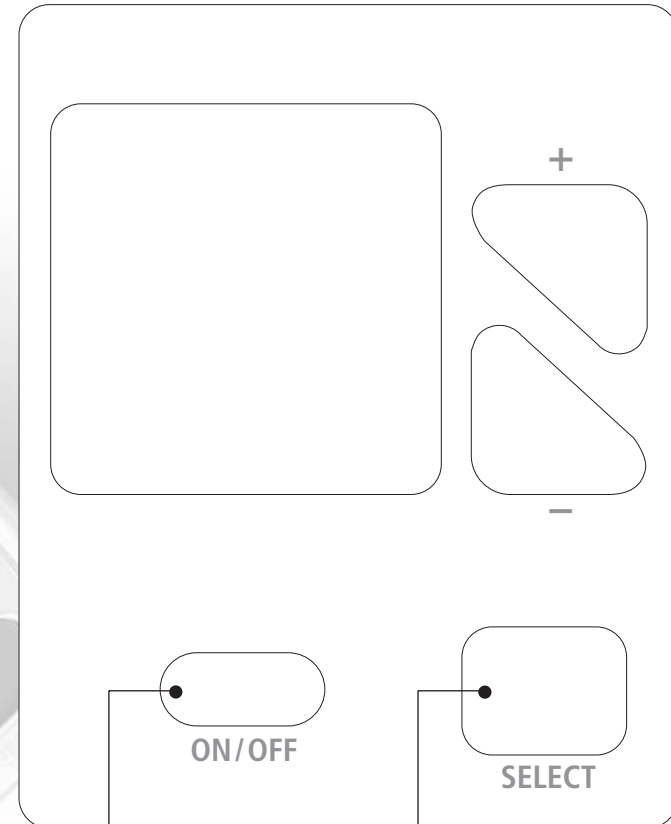ON/OFF:

- Press for 2 sec.
- Until beep sounds

Display:

- Time/Temperature
- Press briefly

Selector switch: Time/Temperature

- Press for 1 sec.
- Select time or temperature
- Select desired value with +/- button
- Lock: press for 2 sec., until the beep sounds (key symbol is displayed)

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